


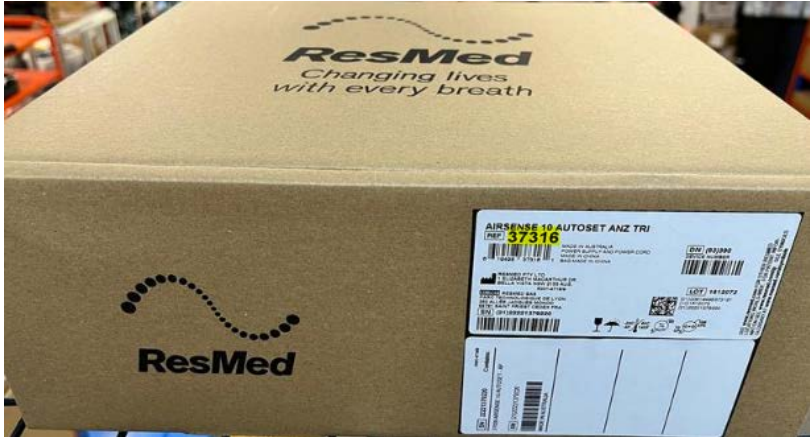


AirSense™ 10 Card-to-Cloud vs CAM Device Comparison

The following information can help assist in physically differentiating between AirSense 10 devices (Elite and AutoSet) that are card-to-cloud only or an AirSense 10 device with a communications access module (CAM) (built-in cellular connectivity).

	AirSense 10 card-to-cloud	AirSense 10 with built-in cellular connectivity
<p>Product codes (visible on the label placed on the outer box packaging)</p>	<ul style="list-style-type: none"> • 37314 – AirSense 10 Elite • 37316 – AirSense 10 AutoSet <p>The product description <i>does not</i> include 3G or 4G in the name.</p> 	<ul style="list-style-type: none"> • 37439 – AirSense 10 Elite 3G • 37436 – AirSense 10 AutoSet 3G • 37460 – AirSense 10 Elite 4G • 37457 – AirSense 10 AutoSet 4G <p>The product description includes 3G or 4G in the name.</p> 

Base codes
(REF no. listed on the back of the device)


Base codes for card-to-cloud devices are different to the device product codes (finished good).

- REF 37021 for device product code 37314
- REF 37029 for device product code 37316



Base codes for 3G & 4G devices are the same as the device product codes (finished good).


- REF 37439 for device product code 37439
- REF 37436 for device product code 37436
- REF 37460 for device product code 37460
- REF 37457 for device product code 37457

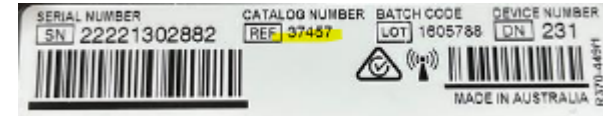

Addendum
(A6 Insert)

378499 – AirSense 10 No-CAM Addendum included with 37316.

In the process of updating 37314 to include the addendum as well for future device builds.

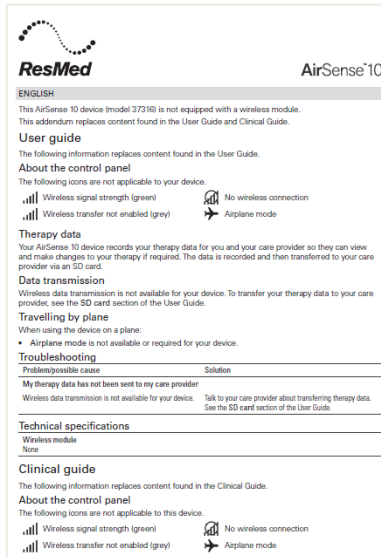


N/A

**myAir™
materials**



No myAir materials included as AirSense 10 card-to-cloud devices are not compatible with myAir.

Note: there was a build of 37314 – AirSense 10 Elite devices where myAir materials were accidentally included.
Serial number range: 22211379569 – 22211383900.

- myAir promo sticker on the top of the device is visible.

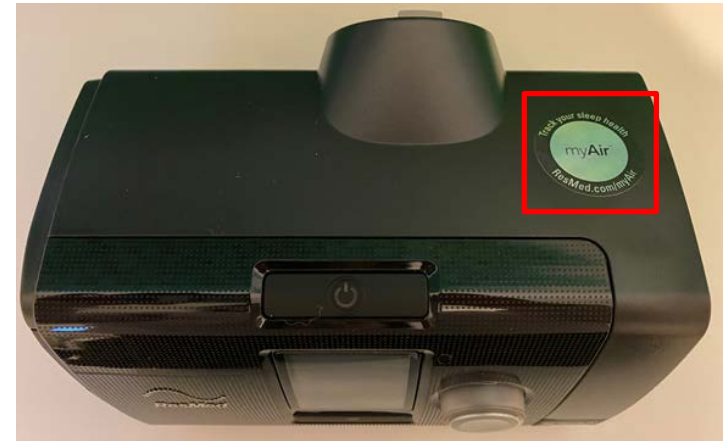

- myAir invite included.

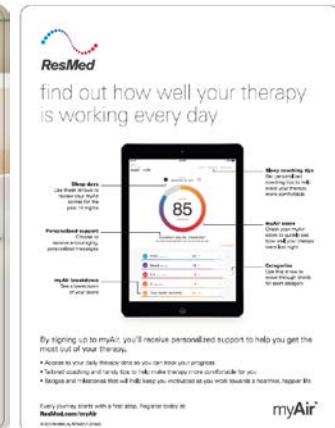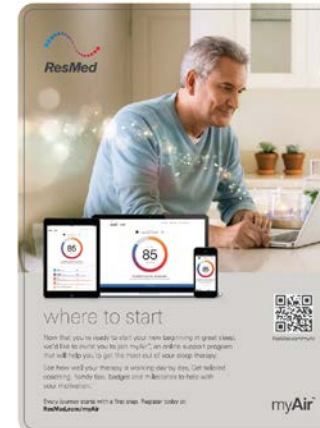


**Non-Ionising
Radiation
Symbol**



Visible, but cannot be used to identify if an AirSense 10 device is card-to-card vs with built-in cellular connectivity.



Visible – required for AirSense 10 devices with built-in cellular connectivity due to the Radiofrequency (RF) radiation.

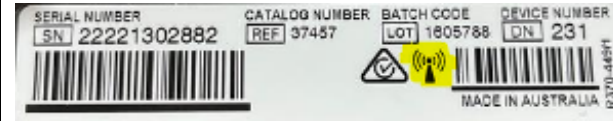

If you have any questions, please contact your ResMed Account Manager.